


222/223

Barrier-Free®
Examination Tables

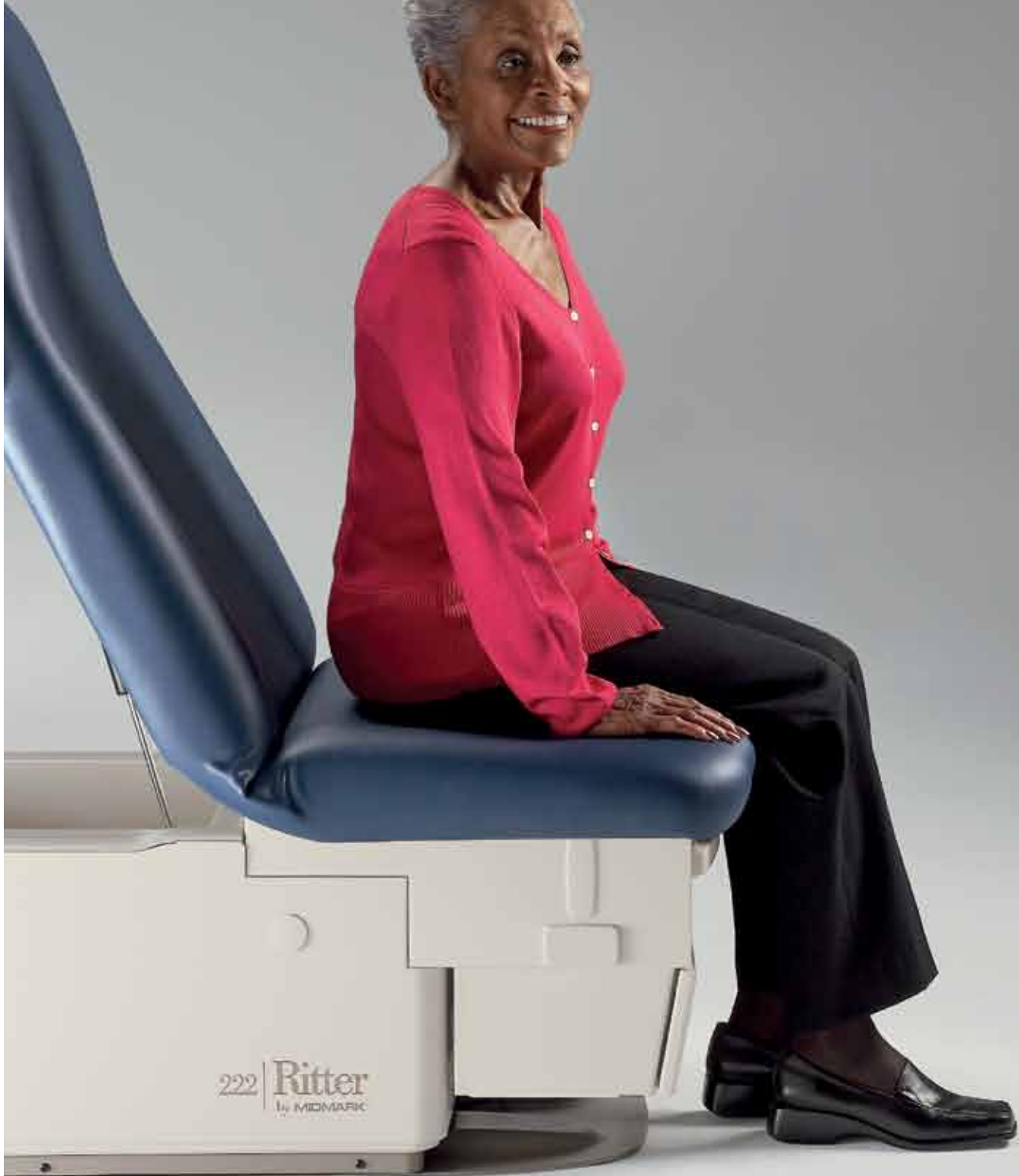

222 | Ritter
by MIDMARK

Ritter®
by MIDMARK

Everyone is entitled to equal access care.

To really care for people, you must understand the patients you care for today...and tomorrow. With an increasingly aging population, your patients may need assistance in accessing the exam table. With an incredible 18-inch low height, Ritter 222 & 223 Barrier-Free® Power Exam Tables make it easy for your practice to accommodate the elderly, expectant mothers or patients with disabilities easily and efficiently. Help them help themselves.



Barrier-Free Exam[®] Tables address the needs of patients and staff...

The Right Thing to Do – For Your Patients and For Your Practice

The days of limited exam rooms for disabled patients are gone. With Ritter Barrier-Free[®] tables, more of your exam rooms can be accessible to your special needs patients. This not only puts your patients at ease, but it also makes your practice more efficient.

Of course, your staff will love the Ritter 222 or 223 more than anyone. An examination table that lowers to 18 inches from the floor means they don't risk injury from lifting patients. Ease of use and accessibility also mean increased efficiency and an overall increase in patient throughput.

Designed in Collaboration with Experts – You

You've been asking us for a truly accessible table. We've taken your advice and done one better: we asked disability experts and disability advocacy groups to offer their opinion, too. Developed in conjunction with these experts, Ritter Barrier-Free[®] tables – the only tables of their kind – reduce the likelihood of distress and injury to patients who might have difficulty moving on and off an industry standard box-style exam table. (The 222 features a manual gas spring assist back and the 223 back section is fully powered.)

The New Industry Standard


Box-style exam tables have been the industry standard in exam rooms for over 50 years; however, they are inaccessible, uncomfortable and inefficient. The choice is clear...Ritter Barrier-Free[®] tables are the new industry standard with powered height and features that are functional and efficient...so, why settle for anything less in your practice? Ritter 222 & 223 Barrier-Free[®] Power Examination Tables are everything you could ask for, at a price you can afford.



The Ritter 222 & 223 allow patients to feel more stable and relaxed while they're waiting or during their consultation as their feet rest comfortably on the floor.

Box tables, even with the step pulled out, suspends patients feet, resulting in an uncomfortable, unstable situation.

222 & 223 Barrier-Free

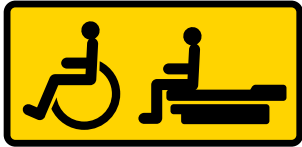

Accessibility and Dignity for All Your Patients

Getting on an exam table can be an anxious, insecure moment, especially for elderly and special needs patients. Ritter 222 & 223 Barrier-Free® Tables lower to 18 inches from the floor, making it easier for patients to get on the table...eliminating excess strain on their backs, as well as, the backs of physicians and staff. This low height means patients can get on the table with no or minimal staff assistance, helping them maintain their dignity.

Measuring 37 inches in the highest position, Ritter Barrier-Free® tables also have the absolute lowest and highest height combination on the market, making exams more comfortable. This maximum low and high height also eliminates excess bending and stooping to help your patients onto the table surface, once again, minimizing stress and strain on the backs of physicians and staff.

Once a patient is on the table, with a touch of the foot control, the Ritter Barrier-Free® tables deliver precise, smooth patient positioning with minimal effort on everyone's part. That goes a long way toward helping your patients – and you – feel relaxed.




The Ritter 222 & 223 were designed with you, your patients and staff in mind. Its 18-inch low height helps the patient move onto and off the table faster, easier and safer.



The 37-inch full height and optional pelvic tilt enables you to conduct an exam or a pap smear without stressing your back.



er-Free® Tables

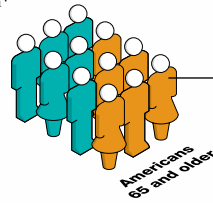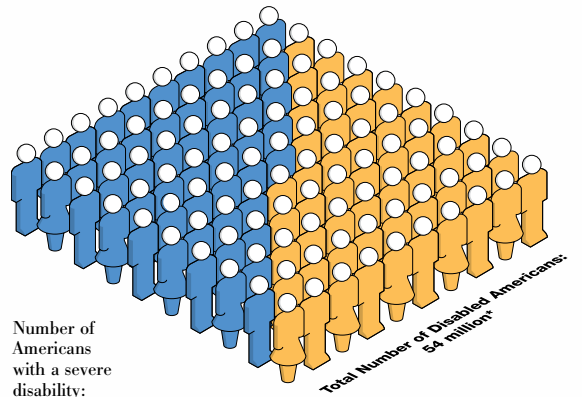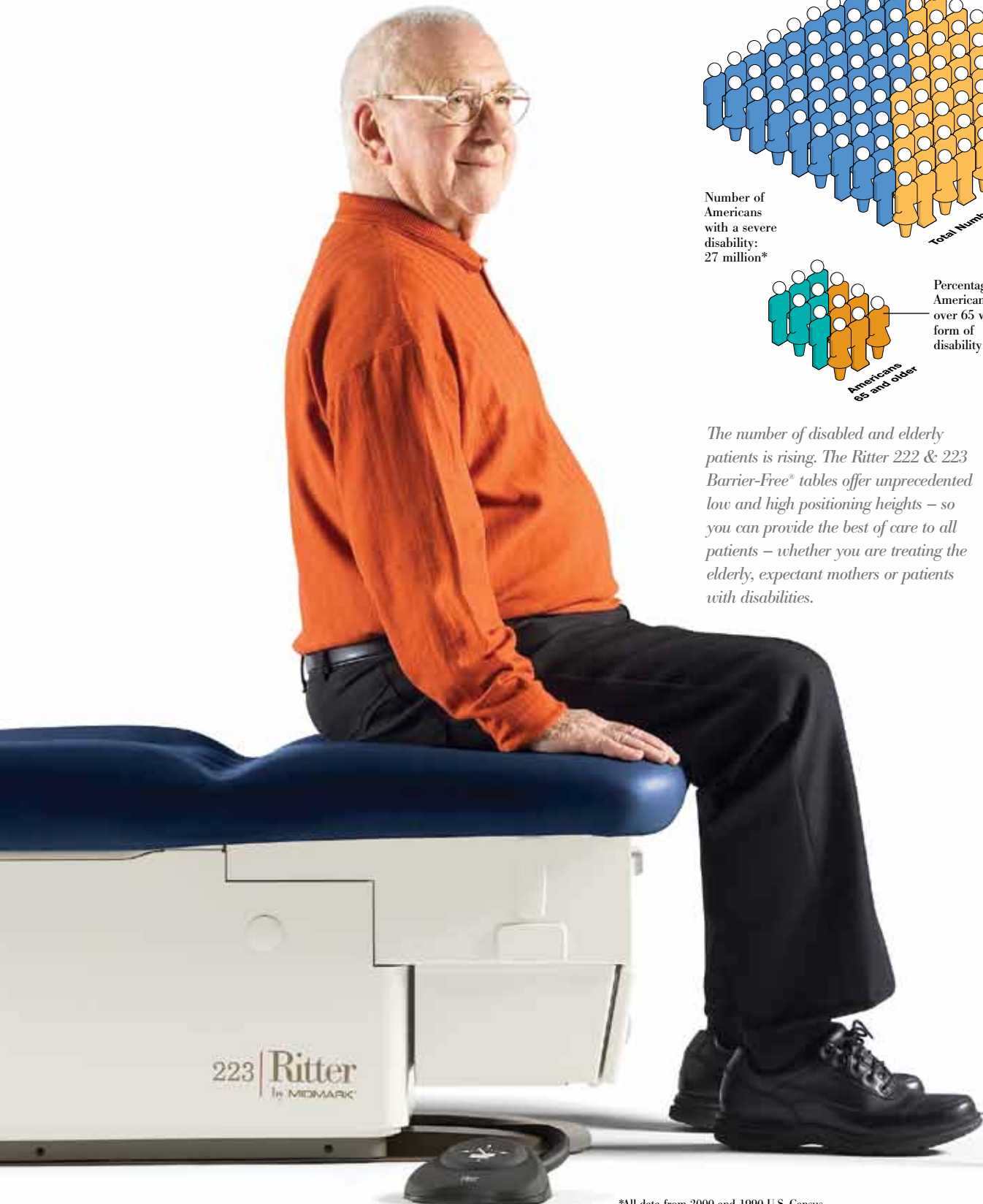

Percentage of Americans over 65 with some form of disability: 50%*

The number of disabled and elderly patients is rising. The Ritter 222 & 223 Barrier-Free® tables offer unprecedented low and high positioning heights – so you can provide the best of care to all patients – whether you are treating the elderly, expectant mothers or patients with disabilities.

*All data from 2000 and 1990 U.S. Census.



Isn't it time to make your exam room more accessible and accommodating for your patients and more efficient for your staff? Midmark's Barrier-Free® examination tables create a safer and more effective environment for patients and medical staff alike.

It's easy to see why Midmark's Barrier-Free® examination tables make the transfer from a wheelchair to the table an easier and safer endeavor.



Both tables offer a stylish, simple to use foot control. The Ritter 223 control allows for simultaneous movement of the base and back section, promoting efficiency.



Removable Upholstery Top
Allows for color change in seconds or additional access for a thorough cleaning. It is easy to remove without the need for tools.

Accessories that add convenience



Articulating Knee Crutches

Unique ball-socket design provides a maximum range of adjustment for optimum knee and leg support.



Flat Upholstery

Provides an excellent working surface for orthopedic, chiropractic or osteopathic uses.



Standard Knee Crutches

Swivel to provide leg support for extended procedures where fatigue is a factor or when the patient lacks sufficient leg strength or control.



Urology Drain Pan

Urology drain pan with hose is ideal for urology procedures requiring fluid and debris removal.



Adjustable Arm System

Offers patients a stable platform in accessing the table and sense of security during positioning. Arms rotate and lock into one of three positions and are easily removable.



Supply Assistant

Provides extra room for exam supplies, such as drapes and gowns, within easy reach.



PC Work Shelf



Supply Storage

Midmark Barrier-Free® tables provide storage for up to four paper rolls. The Supply Assistant rear storage drawer provides ample room for bulky supplies.

Ritter 222 & 223

S P E C I F I C A T I O N S



The Ritter 222 & 223's Exam Assistant™ drawer system provides an area to store basic pelvic and examination supplies at your fingertips. Its unique fold-down front can also turn into a work shelf for exam assistance. This makes the exam easier and more efficient for you.

Ritter 222 Barrier-Free® Power Examination Table

Features

With unprecedented low and high positioning heights, the Ritter 222 has a gas spring assisted back section with an easy access backrest mechanism while the Ritter 223 has a power back system. Both have built-in pillow, hideaway stirrups that lock in one of four lateral positions, foot extension shelf to expand the table length, paper roll holder hidden under headrest, storage for four 21" x 3.5" diameter paper rolls, polystyrene treatment pan, foot control, Exam Assistant™, and leveling screws.

Specifications

Length with footrest extended:
76.7" (194.7 cm)

Height:
Minimum 18" (45.7 cm)
Maximum 37" (94.0 cm)

Patient weight capacity 400#
(181.4 kg)

Paper roll holder (hidden):
Holds 21" x 3.5" paper roll
(53.3 x 8.9 cm)

Electrical receptacle
(115 VAC 5 Amps Max.)
Duplex, hospital grade,
electrical receptacle on
left side

Pelvic tilt (5 degrees)

Drawer heater (35 watt)

Stirrups are standard

Options

Seamless upholstery top
29.7" W x 58.7" L
(75.4 x 149 cm)

Accessories

Rear storage pod (includes
removable divider)
17.8" x 7.0" x 9.8" high
(45.2 x 17.8 x 24.8 cm high)

Adjustable arm system (right
and left side arms included)
Factory installed only

Knee crutches

Procto rest

Welch Allyn® light bracket

Urology pan

Stainless steel pan

Electrical - Ritter 222

115 VAC, 50/60 Hz, 11 Amps
(with heater & receptacle -016)

UL 60601-1 and CAN/CSA C22.2
60601.1-M90

230 VAC, 50/60 Hz, 2.75 Amps
(basic unit -017 excludes
receptacle)

Electrical - Ritter 223

115 VAC, 50/60 Hz, 12.5 Amps
(with heater & receptacle -016)

UL 60601-1 and CAN/CSA C22.2
60601.1-M90



Colors

Available in eight Ritter upholstery colors: Navy, Dusty Blue, Moss, Clay, Pebble Grey, Shadow, Perfect Plum and Black



**CARB 93120.2 Phase II
Compliant**

**Midmark is an ISO 13485 and
ISO 9001 Certified Company.**

For more information
or a demonstration, contact
your Ritter or Midmark
Dealer, or contact Midmark
Corporation, Versailles, Ohio.
Call: 1-800-MIDMARK
Fax: 1-800-365-8631
midmark.com

Ritter®
by MIDMARK

© Midmark Corporation 2002
Products subject to improvement changes without notice.
Litho in U.S.A. 007-0344-00 Rev. G1 (8/11)