HEIGHT ADJUSTABLE BATH TROLLEY

Bolero



Bolero is a height adjustable bath trolley that can be used for transporting, transferring in and out of the bathtub and bathing dependent residents in reclined or fully lying position.

ARJO MOBILITY GALLERY



Designed to support the bathing workflow

The height adjustable stretcher features a padded and thus comfortable lying surface. It enables a single caregiver to safely and efficiently transport a resident directly from the bedside or bathroom into the bath, and back again, without the need for manual lifting or additional transfers. Thus, providing an efficient and more ergonomic alternative to in-bed bathing.

Symmetrical design and foldable mattress enabling flexibility

Bolero has a symmetrical design, enabling either end to be used as a backrest, depending on the orientation of the resident's bed. Along with the foldable mattress facilitating storage, Bolero offers flexibility and ability to be adapted for different resident, caregiver and facility needs.

Supports a range of resident heights and needs

The adjustable backrest enables an upright position to a fully reclined position. The stretcher length can also be shortened by adjusting the head and foot support to optimize the usage for shorter residents.



Key features Bolero comes in two versions, hydraulic or electric.




PRODUCT INFORMATION	BOLERO ELECTRIC	BOLERO HYDRAULIC
Max. patient weight, SWL	136 kg (300 lbs.)	136 kg (300 lbs.)
Total weight of lift (lift + scale)	78 kg (171 lbs.)	69 kg (151 lbs.)
Total width	860 mm (33 7/8")	860 mm (33 7/8")
Total width stretcher	446 mm (17 ½")	446 mm (17 ½")
Total length flat Total length folded up (storage)	1970 mm (77 9/16″) 1000 mm (39 3/8″)	1970 mm (77 9/16″) 1000 mm (39 3/8″)
Min. height of stretcher	515 mm (1/4")	515 mm (1/4")
Max. height of stretcher	1080 mm (42 ½")	1015 mm (40")
Min. height of equipment	1165 mm (45 7/8")	1150 mm (45 ¼″)
Max. height of equipment	1730 mm (68 1/8")	1650 mm (65")
Max. chassis width	1000 mm (39 3/8")	1000 mm (39 3/8")
Min. chassis width between legs	499 mm (19 5/8")	499 mm (19 5/8")
Max. chassis length	823 mm (32 13/32")	823 mm (32 13/32")
Turning radius (flat)	2030 mm (79 7/8")	2030 mm (79 7/8")
Lift stroke range	565 mm (22 ¼″)	500 mm (19 5/8″)

ACCESSORIES		
Safety belts	To ensure resident safety. Bolero comes with two safety belts – one chest belt and one hip belt. Use according to resident need.	
Batteries and charger 24V	Quick-replaceable battery system with two batteries, ensuring availability for use at all times.	
Side support (optional)	To reduce risk of resident falling out of stretcher, and for resident comfort.	
Scale (optional)	Foot-operated scale to weigh resident.	

- Blue area shows the minimum working area required for the staff to be able to use the mechanical aids in an ergonomic way from one side.
- Light blue area shows required extension of working area to facilitate activities from either side to provide adequate access.

For more info on installation requirements visit Architect portal at www.arjo.com

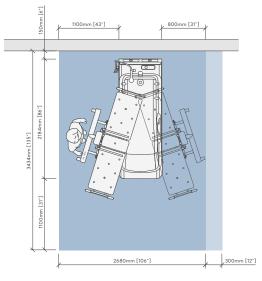

ario

March 2023. Only Arjo designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by Arjo. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice. ® and ™ are trademarks belonging to the Arjo group of companies. Arjo AB • Hans Michelsensgatan 10 • 211 20 Malmö • Sweden • +46 10 335 4500