# **Concerto and Basic**



Concerto and Basic shower trolley help to facilitate efficient showering and hygiene care routines for dependent patients or residents who need to remain in a lying position.

Promotes correct working posture.

#### **ARJO MOBILITY GALLERY**





## Promotes correct working posture

Both Concerto and Basic are height-adjustable, allowing the caregiver to reach their preferred working height. The height-adjustability is operated with a hydraulic pump. Concerto is also available as an electrically powered model.

#### **Enables comfort and support**

To enhance comfort, provide support and aid in positioning the resident, Concerto and Basic come with a padded mattress, head cushion and wedge pillow. Concerto also offers a back support which is adjustable to three positions. Concerto can also be tilted sideways for convenient cleaning and drying.

### Ease of Manouvreability

Concerto offers straight steering on two castors for transporting along corridors, and automatic horizontal levelling when lowered against a bed.



# **Key features**

#### SIDE SUPPORTS


Lockable and foldable side supports that can be lowered to facilitate lateral transfers

**HIGHT ADJUSTABLE**Features pedal and hydraulic

actuator, for easy height

adjustability

pump or electrically powered

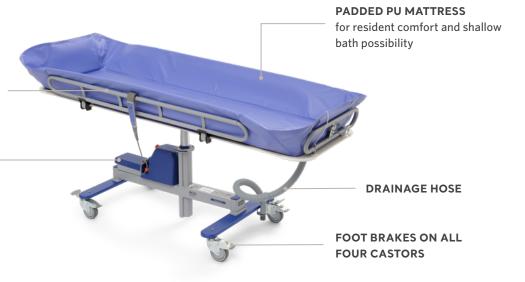

Max total safe working load 175 kg (385 lbs) 150 kg (330 lbs) Max resident/patient weight 320 mm (12 1/2 in) hydraulic Lifting stroke 300 mm (11 3/4 in) electric 740-1060 mm (29-41 3/4 in) hydraulic Height, top position of stretcher 740-1040 mm (29-41 in) electric 72 kg (159 lbs) Basic Weight 75 kg (165 lbs) Concerto 79 kg (174 lbs) electric Width, stretcher 620 mm (23 13/32 in) Total width (including side supports) 755 mm (29 3/4 in) Total length (different stretcher lengths) 1600 mm, 1850 to 1950 mm, 2250 mm (63 in, 73 to 76 3/4 in, 88 1/2 in) Length of chassis (including wheel frame) 1155 mm (45 1/2 in) Height of chassis 160 mm (6 1/4 in) +10 to +40 C (+50 to +104 F) Operating Temperature -20 to +70 C (-4 to +158 F) Transport -20 to +70 C (-4 to +158 F) Storage 30-75% Operating

10-80% Transport

10-80% Storage

### **SPACE REQUIREMENTS**

Humidity

- Blue area shows the minimum working area required for the staff to be able to use the mechanical aids in an ergonomic way from one side.
- Light blue area shows required extension of working area to facilitate activities from either side to provide adequate access.

For more info on installation requirements visit Architect portal at www.arjo.com

