


PRESSURE INJURY PREVENTION

Simulflex™


Foam mattress range

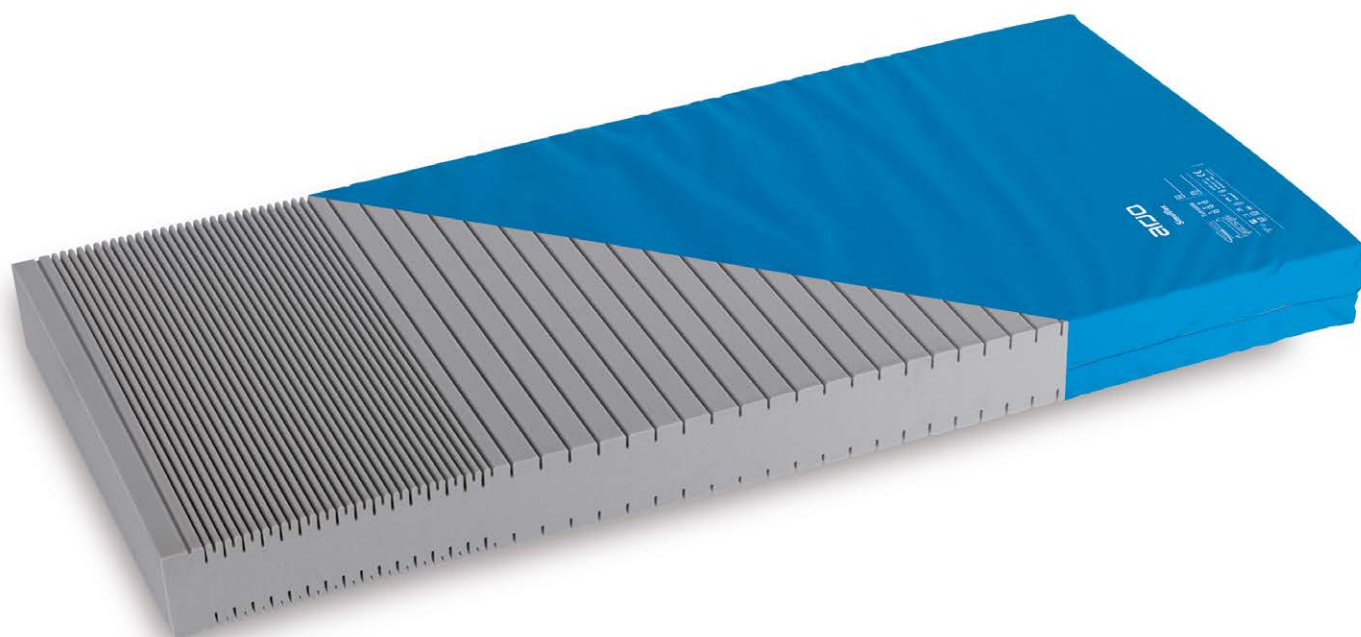

The Simulflex range has been developed to provide a 24 hour care solution. The range includes a mattress replacement, mattress overlay and seat cushion.


Simulflex™ Range

This range is made of high quality polyurethane foam and available in variable sizes. In addition, the overlay includes corner securing straps for extra security.

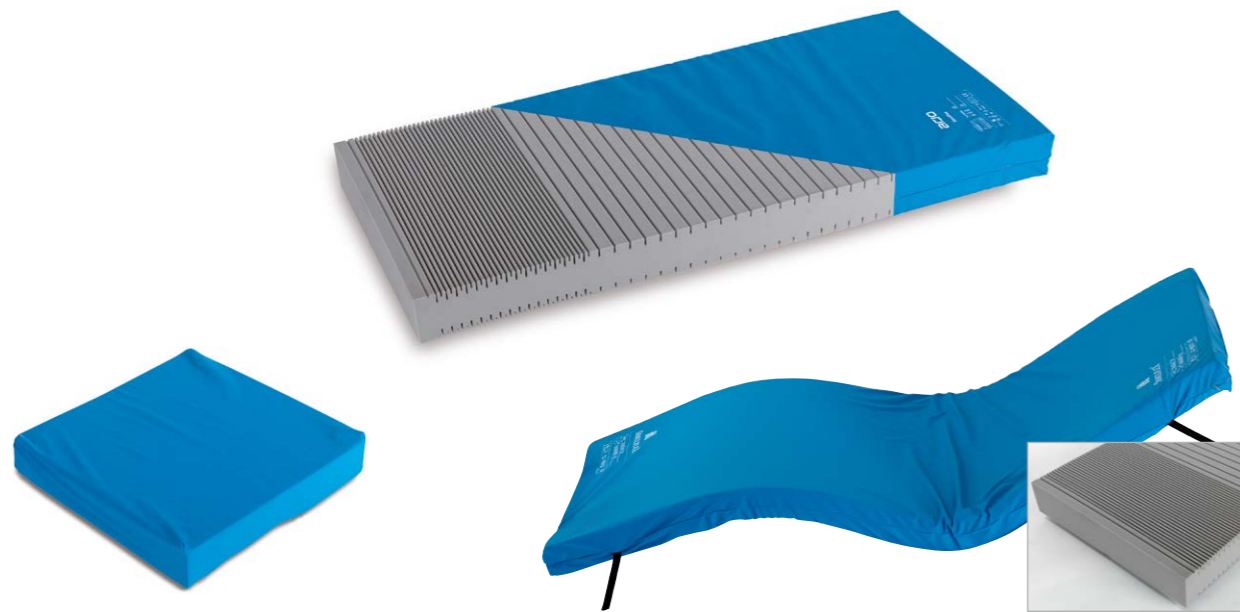


24 Hour Care Package - Seat cushion available



Overlay System - A mattress overlay is also available and is easy to install

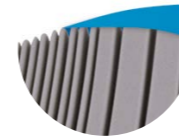


ANATOMICALLY INSPIRED ZONES
Features horizontally cut foam to provide 3 zones of pressure redistribution.

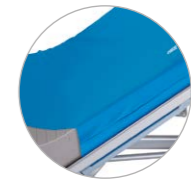

CONCEALED ZIP
The cover has a concealed zip designed to help improve infection prevention protocols.

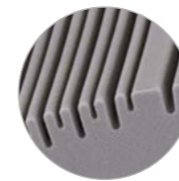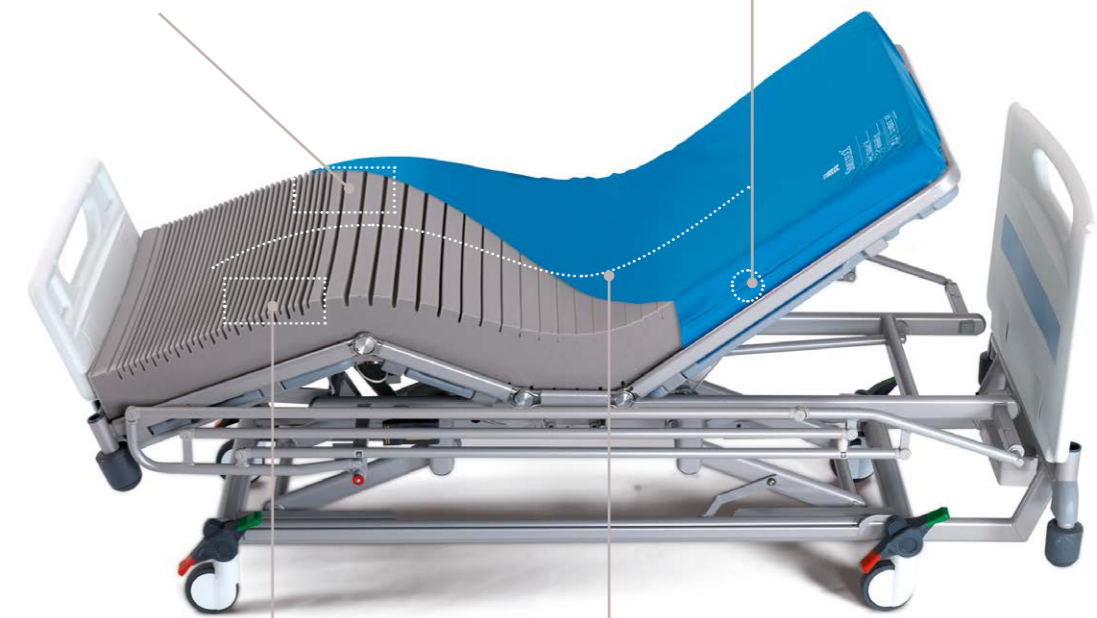

HEAD AND HEEL ZONE
The head and heel zones include smaller cuts in the foam designed to improve comfort for the patient.

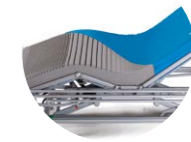

SUITABLE FOR PROFILING BEDS
Suitable for use in the gatched or supine position.

Comfort and safety in mind


Suitable for use in a variety of care environments including hospitals, and community care

Safe and convenient for patients and carers

Simulflex range has been designed to provide reliable pressure redistribution performance.

Cover features

- Two way stretch
- Vapour permeable
- Water resistant
- Infection Control is aided by Antimicrobial properties in the cover material
- Easy to clean and maintain
- Concealed zip
- Clearly marked covers with TurnTable™ guide and washing instructions

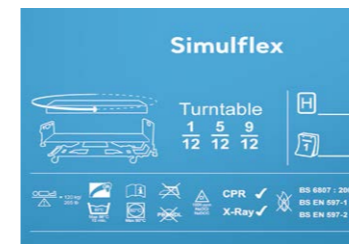

TurnTable Guide - All Simulflex mattress replacement systems carry the TurnTable Guide

Specifications*

4-WAY TURN MATTRESS REPLACEMENT	DIMENSIONS H X W X D**	ORDER CODE
MR Standard	1980 x 860 x 150 mm (78 x 33.9 x 5.9")	120SX198/086
MR Community	1900 x 920 x 150 mm (74.8 x 36.2 x 5.9")	120SX190/092
MR Contoura	1820 x 780 x 150 mm (71.7 x 30.7 x 5.9")	120SX182/078
Max. patient weight	120 kg (264 lbs)	
Lifetime expectancy: 5 years for foam and 2 years for cover		
Overlay system		
OL Standard	1980 x 860 x 85 mm (78 x 33.9 x 33.5")	SO198/086
OL Community	1900 x 920 x 85 mm (74.8 x 36.2 x 33.5")	SO190/092
Seat cushion		
Seat	430 x 430 x 70 mm (6.9 x 6.9 x 2.8")	SS043/043
User weight limit:	90 kg (198 lbs)	
Mattress replacement fire retardancy		
MR System meets BS EN 597 -1: 1995, BS EN 597 -2: 1995 and BS 6807: 1996, Section 2, Ignition Sources 0,1 and 5 Overlay System meets BS EN 597-1: 1995, BS EN 597 -2: 1995 and BS 6807: 1996, Section 2, Ignition Sources 0,1 and 5		
X-Ray	Tested to FDA 21 CFR Ch. 1 (4-1-95 Edition) HHS & CEI/IEC 601-1-3:1994-07 (X-ray attenuation)	
CPR	Chest compression for CPR purposes is not impeded	

*Specifications are subject to change without notice

**Simulflex mattresses are available in other custom made sizes. All sizes must be specified in increments of 2 cm

October 2019 Only Arjo designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by Arjo. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice. ® and ™ are trademarks belonging to the Arjo group of companies.
© Arjo, 2019

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure injury and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind.

Arjo AB · Hans Michelsensgatan 10 · 211 20 Malmö · Sweden · +46 10 335 4500

www.arjo.com